

TLS 01

Thoraco lumbo sacral support

Ref. SRT 117

Indications

It is used in case of intervertebral disc herniation, back pain, postoperative, osteoporosis, for stabilization of the thoracolumbosacral spine.

Description

The support is made of elastic fabric, provided with:

- 2 lateral plastic splints for support
- 2 long, paravertebral plastic splints for support
- 2 elastic straps that offer an increased compression on the lumbar region
- 2 elastic straps that surround the shoulders and are fixed on the chest.

The straps are provided with sleeves for armpit protection.

Hook and loop fastener on the front side.

Width: front 16 cm, back height between 48 - 53 cm, depending on size.


For an easy application it is provided with helping strips.

The 2 lateral and the 2 paravertebral plastic rods stabilize the thoracolumbar region, slightly reducing lordosis and kyphosis, and the compression exerted by the elastic bands stimulates local metabolism, reducing muscle tension and pain.


The 2 elastic bands that surround the shoulders and are fixed on the chest, correct by traction the position of the trunk. It is a traditional elastic corset that, due to the material, allows normal thoraco abdominal breathing, while ensuring uniform compression on the entire covered surface.

Composition

Elastic fabric: 50% polyester, 33% polyamide, 17% elastane, with a structure that allows high air permeability in both directions.

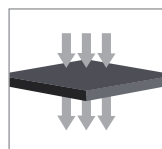

air transfer

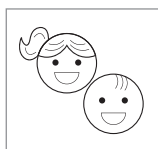

also available in Junior version




- Available colors: black, beige



Product details



- hook and loop fastening on the front
- sleeves for axillary protection

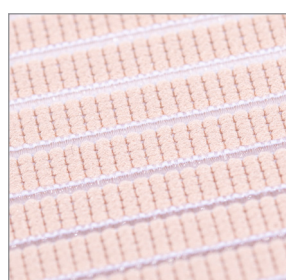

- elastic fabric with a structure that allows high air permeability in both directions

Available sizes

| Sizes | 1 | 2 | 3 | 4 | 5 | 6 |
|------------------------------------|----------|----------|----------|----------|----------|----------|
| Waist circumference (cm) | 65-75 | 75-85 | 85-95 | 95-105 | 105-115 | 115-125 |
| Height of support on the back (cm) | 48 | 49 | 50 | 51 | 52 | 53 |
| Reference | SRT117M1 | SRT117M2 | SRT117M3 | SRT117M4 | SRT117M5 | SRT117M6 |

The circumference of the waist is measured

