


MALLEOSPORT

Mobile ankle brace with elastic strap

Ref. SRT 347

INDICATIONS

It is used in cricks, luxations, post fractures, inflammatory processes, after the removal of the orthopedic plaster cast.

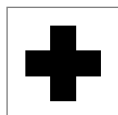

MEDICAL



SPORT



DAILY

DESCRIPTION

Orthosis with anatomical design, made of textile material, provided with:

- 2 straps with hook and loop fastener
 - one elastic strap for the adjustment of the compression.
- The heel area is open.


Hook and loop fastener on the front side.

Height: 21 cm.


The orthosis is an ankle support, ensuring stability while wearing.

COMPOSITION

Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

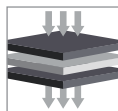

AIR FLOW

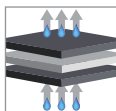

WATER
TRANSFER

AVAILABLE COLOURS

Black



AVAILABLE SIZES

Universal size

For ankle circumference between 20 – 30 cm.

