

GENU REHA PO

Mobile knee brace with adjustable articulations

Ref. SRT 309

Indications

It is used for postoperative recovery, control of the flexion and extension of the knee joint.

Description

The orthosis is open in the front, made of textile, provided with:

- 2 metal splints with mobile and adjustable articulations, one on each side. The adjustment device is easily accessed on the sides, covered with a plastic cap. The angle is adjusted with 2 rivets. The possibilities of limiting the movement are:

- **flexion** 15°, 30°, 50°, 70°, 90°, 110°

- **extension** 15°, 30°, 50°, 70°, 90°, 110°

- 2 hook and loop straps on the leg


- 2 hook and loop straps on the thigh.

Length: 44 cm.


The orthosis is a support for the knee joint, the 2 lateral metal splints providing stability and safety in use. Due to its particular shape, it can be used also in case of a deformed knee. The compression can be adjusted by the frontal closure system and the 4 fastening straps.

Composition

Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

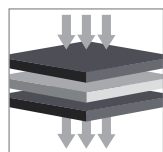

air transfer

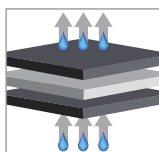

water transfer



Product details



- the orthosis is open in the front
- hook and loop straps



- the popliteal area is free




- metal splints with mobile, adjustable joint



- Available color: black

Available sizes

Sizes	1	2	3	4	5	6	7	8
Knee circumference (cm)	33-35	35-37	37-39	39-41	41-43	43-45	45-47	47-49
Reference	SRT309M1	SRT309M2	SRT309M3	SRT309M4	SRT309M5	SRT309M6	SRT309M7	SRT309M8

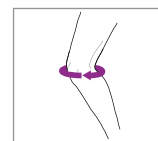

The circumference of the knee is measured at the middle of the knee cap