

TRIAGEN EXTRA

Mobile knee brace with patellar support and lateral articulations

Ref. SRT 341

Indications

It is used in lateral instability of the knee, traumas or injuries at the level of the knee, inflammatory processes, for postoperative recovery.

Description

The orthosis is made of textile, open in the front, provided with:

- a non-slip silicone elastic band (10mm) inside the upper part
- 2 metal splints with mobile articulations, one on each side
- 2 hook and loop straps on the leg
- 2 hook and loop straps on the thigh
- cutout in the popliteal region for comfort.

At the level of the opening for the knee cap, the textiles are overlapped, and after the orthosis is closed, a circular patellar support is created.


Length: 38 cm.

Due to its particular shape, it can be used also in case of a deformed knee.

The compression can be adjusted by the frontal closure system and the 4 fastening straps.

Composition

Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

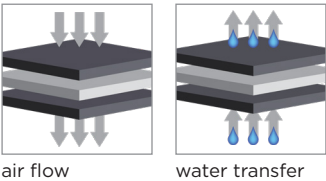

air flow

water transfer



- Available colour: black

Available sizes

Sisez	1	2	3	4	5	6	7	8
Knee circumference (cm)	33-35	35-37	37-39	39-41	41-43	43-45	45-47	47-49
Reference	SRT341M1	SRT341M2	SRT341M3	SRT341M4	SRT341M5	SRT341M6	SRT341M7	SRT341M8

The circumference of the knee is measured at the middle of the knee cap



Product details



- orthosis open in the front
- hook and loop fastener
- non-slip silicone elastic band



- the popliteal area is free



- metal splints with mobile joint

