


TRIAMED

Shoulder elbow wrist hand brace

Ref. SRT 219

INDICATIONS


It is used in case of humerus sprain, humerus collar fracture, the humeral tuberosity fracture without movement, shoulder blade fracture, scapulohumeral dislocation.

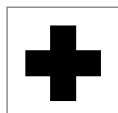

MEDICAL



SPORT

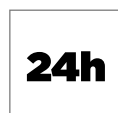

DAILY

DESCRIPTION

Orthosis with anatomical design, made of textile, provided with 2 components:

- one arm support with 2 straps adjustable by cutting them and repositioning the "crocodile" type hook and loop fastening


- one adjustable thoracic belt to be fixed around the thorax and arm after the immobilization of the arm. The regulation is done by cutting at the level of the stitches to avoid fabric breaking.

Hook and loop fastener.


The orthosis is a support for the shoulder, arm and the forearm. Due to the thoracic belt the support does not allow the abduction of the arm and external rotation.

COMPOZIȚIE

Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture outward.

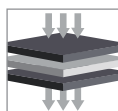

AIR FLOW

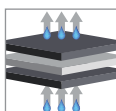

WATER
TRANSFER

AVAILABLE COLOURS

Black



AVAILABLE SIZES

Universal size

For arm and thorax circumference between 75 – 147 cm.

