

## TRIAGEN PATELLA

Mobile knee brace

Ref. SRT 342

### Indications

It is used in cricks, luxations, post consolidated fractures, inflammatory processes, lateral knee joint instability, instability of the knee cap, osteoarthritis. It is recommended both as orthopedic treatment and postoperative.

### Description

It is a tubular, closed orthosis, made of textile, provided with:

- a non-slip silicone elastic band (10mm) inside the upper part
- 4 spring type metal splints, 2 on each side
- 2 pieces for the patella, not attached to the orthosis, in the shape of the letter "Y", one for use on the left limb, the other for the right limb. The patella pieces are made of neoprene and are provided with hook and loop fastener on all 3 edges
- patellar support (polyethylene foam).

Length: 29 cm.

The orthosis is a support for the knee joint, the 2 spring type metal splints for each side providing stability and safety in use. Due to its tubular, sleeve shape, it is not recommended to be used in case of oedema or significant swelling.

### Composition

Neopren coated on both sides with fabric, 100% polyester.




### Product details



- tubular, closed orthosis
- 2 neopren pieces provided with hook and loop fastener on all 3 edges, one for the right leg and one for the left leg



- patellar support made of polyethylene foam

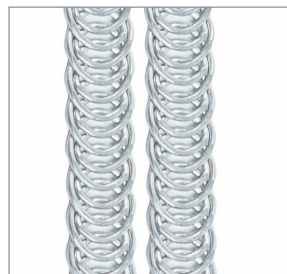

- spring type metal splints



- Available colour: black

### Available sizes

Size	1	2	3	4	5	6	7	8
Knee circumference (cm)	33-35	35-37	37-39	39-41	41-43	43-45	45-47	47-49
Reference	SRT342M1	SRT342M2	SRT342M3	SRT342M4	SRT342M5	SRT342M6	SRT342M7	SRT342M8

The circumference of the knee is measured at the middle of the knee cap

