

## TRIACROSS

Mobile knee brace

Ref. SRT 317

### Indications

It is used in cricks, luxations, post consolidated fractures, lateral knee joint instability, inflammatory processes, anterior cruciate ligament injury. It is recommended both as orthopedic treatment and postoperative.

### Description

It is a tubular, closed orthosis, made of textile, provided with:

- a non-slip silicone elastic band (10mm) inside the upper part
- 4 spring type metal splints, 2 on each side
- patellar support (polyethylene foam)
- elastic strap crossed in the front, under the patellar support, and on the back, in the popliteal region.

Length: 29 cm.

The orthosis is a support for the knee joint, the 2 spring type metal splints for each side providing stability and safety in use.

Due to its tubular, sleeve shape, it is not recommended to be used in case of oedema or significant swelling.

### Composition

Neopren coated on both sides with fabric, 100% polyester.




### Product details



- tubular, closed orthosis



- elastic band crossed in the front, under the patellar support

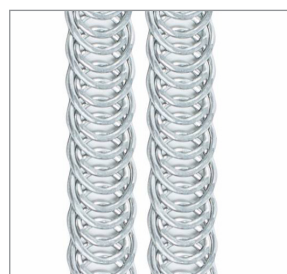

- spring type metal splints



- Available colour: black

### Available sizes

Sizes	1	2	3	4	5	6	7	8
Knee circumference (cm)	33-35	35-37	37-39	39-41	41-43	43-45	45-47	47-49
Reference	SRT317M1	SRT317M2	SRT317M3	SRT317M4	SRT317M5	SRT317M6	SRT317M7	SRT317M8

The circumference of the knee is measured at the middle of the knee cap

