

SENIOR 01

Back support

Ref. SRT 120

Indications

It is used in case of incipient intervertebral disc herniation, moderate low back pain, postoperative, in case of osteoporosis, for moderate stabilization of the lumbar spine.

Description

The support is made of textile, provided with 4 integrated plastic rods for support.

Hook and loop fastener on the front side.


Width: front 16 cm, back 27 cm.

For an easy application it is provided with helping strips.


The 4 plastic rods stabilize the lumbar region, slightly reducing lordosis, and the compression exerted by the elastic bands stimulates local metabolism, reducing muscle tension and pain.

Composition

Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

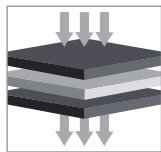

air transfer

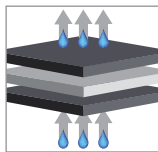

water transfer



Product details



- helping strips for an easy application




- hook and loop fastener on the front side



- Available color: black

Available sizes

Sizes	1	2	3	4	5	6
Waist circumference (cm)	65-75	75-85	85-95	95-105	105-115	115-125
Reference	SRT120M1	SRT120M2	SRT120M3	SRT120M4	SRT120M5	SRT120M6

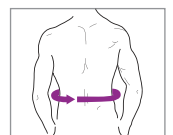

The circumference of the waist is measured