

GEROMED 02

Mobile knee brace

Ref. SRT 340

Indications

It is used in cricks, luxations, post consolidated fractures, inflammatory processes, lateral knee joint instability. It is recommended both as orthopedic treatment and postoperative.

Description

The orthosis is open at the back, made of textile, provided with:

- a non-slip silicone elastic band (10mm) inside the upper part
- 4 spring type metal splints, 2 on each side
- patellar support (polyethylene foam)
- 3 neoprene fastening straps: 2 front and one rear closure
- 2 straps of elastic fastening material. The elastic fastening follows the conical shape of the thigh, preventing the orthosis from slipping.

Length: 32 cm.

The orthosis is a support for the knee joint, the 4 spring type metal splints, 2 for each side providing stability and safety in use. Due to its particular shape, it can be used also in case of a deformed knee. The compression can be adjusted by the back closure system.

Composition

Neopren coated on both sides with fabric, 100% polyester.




Product details



- orthosis open on the back side
- hook and loop fastener



- patellar support made of polyethylene foam

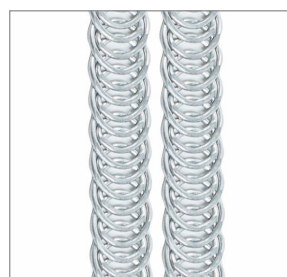

- spring type metal splints



- Available colour: black

Available sizes

Sizes	1	2	3	4	5
Knee circumference (cm)	33-36	36-39	39-42	42-45	45-49
Reference	SRT340M1	SRT340M2	SRT340M3	SRT340M4	SRT340M5

The circumference of the knee is measured at the middle of the knee cap

